



A Sysmex Group Company



PML/RAR α (RARA) Translocation, Dual Fusion Probe


REF: LPH 023-A / LPH 023-A50

Analyte Specific Reagent: Analytical and performance characteristics are not established.

Fluorescence *in situ* hybridisation (FISH) is a technique that allows the visualisation of DNA sequences upon chromosomes. The technique uses DNA probes that hybridise to entire chromosomes or single unique sequences, and serves as a powerful adjunct to classic cytogenetics. Recent developments have meant that this valuable technique can now be applied as an essential tool in prenatal, haematological and pathological chromosomal analysis. Target DNA, after fixation and denaturation, is available for annealing to a similarly denatured, fluorescently labelled DNA probe, which has a complementary sequence. Following hybridisation, unbound and non-specifically bound DNA probe is removed and the DNA is counterstained for visualisation. Fluorescence microscopy then allows the visualisation of the hybridised probe on the target material.

Probe Specification

PML, 15q24.1, Red
RAR α , 17q21.1-q21.2, Green

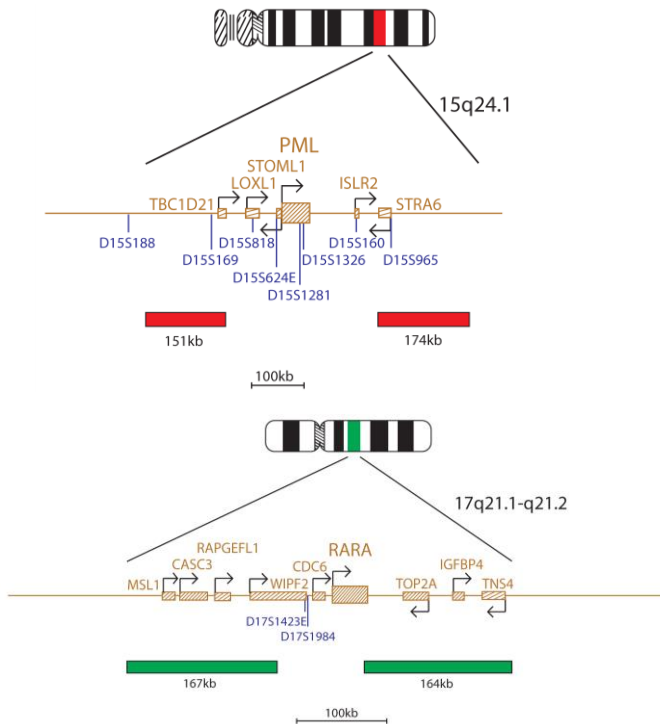

The PML probe mix, labelled in red, consists of a 151kb probe centromeric to the PML gene and a 174kb probe telomeric to the PML gene. The RAR α (RAR α) probe mix, labelled in green, consists of a 167kb probe centromeric to the RAR α gene, including the CASC3 gene, and a 164kb probe, including the telomeric end of the RAR α gene as well as the TOP2A and IGFBP4 genes.

Materials Provided

Probe: 100 μ l per vial or 500 μ l per vial
Probe concentration: Amount of red probe 1.40-2.36ng/ μ l
Amount of green probe 5.48-8.20ng/ μ l

The probe is provided in hybridisation solution (Formamide; Dextran Sulphate; SSC) and is ready to use.

Warnings and Precautions

1. For professional use only.
2. Wear gloves when handling DNA probes.

3. Probe contains formamide, which is a teratogen; do not breathe the fumes or allow skin contact. Wear gloves, a lab coat, and handle in a fume hood. Upon disposal, flush with a large volume of water.
4. Dispose of all hazardous materials according to your institution's guidelines for hazardous waste disposal.
5. Operators must be capable of visually distinguishing between red, blue and green.

Storage and Handling

The kit should be stored between -25°C to -15°C in a freezer until the expiry date indicated on the kit label. Wear gloves, a lab coat, and handle in a fume hood. Ensure that exposure of the probe to laboratory lights is limited at all times.

Known Cross-Reactivity

No known cross-reactivity.

Additional Information

For additional product information please contact the CytoCell Technical Support Department.

T: +44 (0)1223 294048
E: techsupport@cytoCELL.com
W: www.ogt.com

Labelling according to Regulation (EC) No. 1272/2008 [CLP]

Hazard pictograms (CLP):



GHS07

GHS08

Signal word (CLP): Danger

Hazardous ingredients: Formamide <100%

Hazard statements (CLP):

H315 – Causes skin irritation
H319 – Causes serious eye irritation
H360 – May damage fertility or the unborn child

Precautionary statements (CLP):

P202 – Do not handle until all safety precautions have been read and understood
P280 – Wear eye protection, protective clothing, protective gloves
P302+P352 – IF ON SKIN: Wash with plenty of soap and water
P305+P351+P338 – IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing
P308+P313 – IF exposed or concerned: Get medical advice/attention
P362+P364 – Take off contaminated clothing and wash it before reuse
P501 – Dispose of contents/container to hazardous or special waste collection point, in accordance with local, regional, national and/or international regulation

Refer to Safety Data Sheet for more information.

Patents and Trademarks

CytoCell is a registered trademark of CytoCell Ltd.
This product contains technology licensed from Life Technologies Corporation that is available for human diagnostics or life science research use only.

CytoCell Ltd.

Oxford Gene Technology,
418 Cambridge Science Park,
Milton Road,
Cambridge, CB4 0PZ, UK
T: +44(0)1223 294048
F: +44(0)1223 294986
E: probes@cytoCELL.com
W: www.ogt.com

