



A Sysmex Group Company



## TEL/AML1 (ETV6/RUNX1) Translocation, Dual Fusion Probe


REF: LPH 012-A / LPH 012-A50

**Analyte Specific Reagent: Analytical and performance characteristics are not established.**

Fluorescence *in situ* hybridisation (FISH) is a technique that allows the visualisation of DNA sequences upon chromosomes. The technique uses DNA probes that hybridise to entire chromosomes or single unique sequences, and serves as a powerful adjunct to classic cytogenetics. Recent developments have meant that this valuable technique can now be applied as an essential tool in prenatal, haematological and pathological chromosomal analysis. Target DNA, after fixation and denaturation, is available for annealing to a similarly denatured, fluorescently labelled DNA probe, which has a complementary sequence. Following hybridisation, unbound and non-specifically bound DNA probe is removed and the DNA is counterstained for visualisation. Fluorescence microscopy then allows the visualisation of the hybridised probe on the target material.

### Probe Specification

TEL (ETV6), 12p13.2, Red  
AML1 (RUNX1), 21q22.1, Green

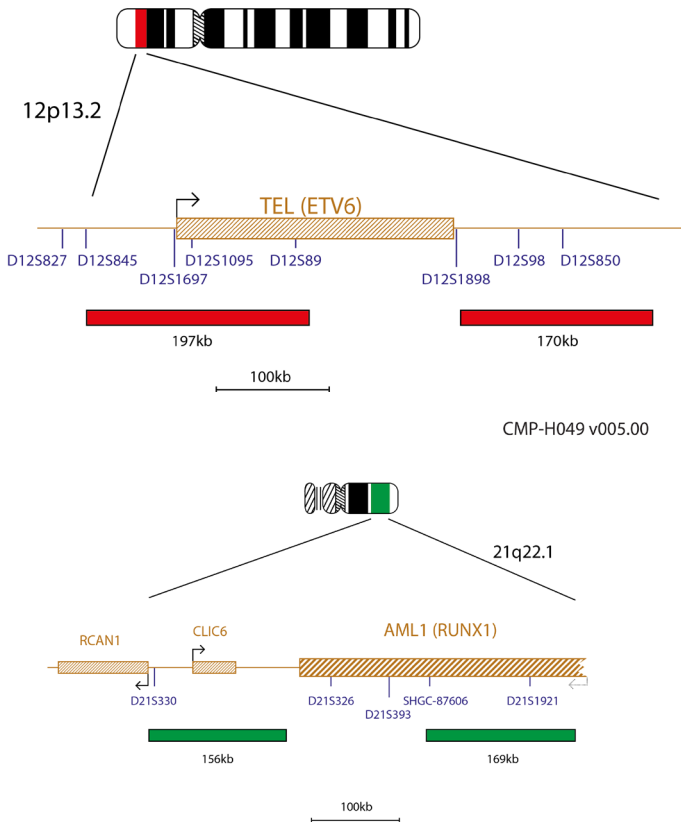

The TEL probe mix, labelled in red, contains a probe covering a 197kb region at the telomeric end of TEL (ETV6) and a second probe that extends 170kb in a centromeric direction beyond the TEL (ETV6) gene. The AML1 (RUNX1) probe set consists of two probes labelled in green. One locates to a 156kb region centromeric to the AML1 (RUNX1) gene and covers the CLIC6 gene. The other probe covers a 169kb region in the AML1 (RUNX1) gene that includes the markers SHGC-87606 and D21S1921.

### Materials Provided

**Probe:** 100µl per vial or 500µl per vial

**Probe concentration:** Amount of red probe 4.55-7.67ng/µl

Amount of green probe 17.8-26.7ng/µl

The probe is provided in hybridisation solution (Formamide; Dextran Sulphate; SSC) and is ready to use.

### Warnings and Precautions

1. For professional use only.
2. Wear gloves when handling DNA probes.
3. Probe contains formamide, which is a teratogen; do not breathe fumes or allow skin contact. Wear gloves, a lab coat, and handle in a fume hood. Upon disposal, flush with a large volume of water.
4. Dispose of all hazardous materials according to your institution's guidelines for hazardous waste disposal.
5. Operators must be capable of visually distinguishing between red, blue and green.

### Storage and Handling

The kit should be stored between -25°C to -15°C in a freezer until the expiry date indicated on the kit label. Store the probe vial in the dark. Ensure that exposure of the probe to laboratory lights is limited at all times.

### Known Cross-Reactivity

No known cross-reactivity.

### Additional Information

For additional product information please contact the CytoCell Technical Support Department.

T: +44 (0)1223 294048

E: techsupport@cytozell.com

W: www.ogt.com

### Labelling according to Regulation (EC) No. 1272/2008 [CLP]

Hazard pictograms (CLP):



GHS07

GHS08

Signal word (CLP): Danger

Hazardous ingredients: Formamide <100%

Hazard statements (CLP):

H315 – Causes skin irritation

H319 – Causes serious eye irritation

H360 – May damage fertility or the unborn child

Precautionary statements (CLP):

P202 – Do not handle until all safety precautions have been read and understood

P280 – Wear eye protection, protective clothing, protective gloves

P302+P352 – IF ON SKIN: Wash with plenty of soap and water

P305+P351+P338 – IF IN EYES: Rinse cautiously with water for several minutes.

Remove contact lenses, if present and easy to do. Continue rinsing

P308+P313 – IF exposed or concerned: Get medical advice/attention

P362+P364 – Take off contaminated clothing and wash it before reuse

P501 – Dispose of contents/container to hazardous or special waste collection point, in accordance with local, regional, national and/or international regulation

Refer to Safety Data Sheet for more information.

### Patents and Trademarks

CytoCell is a registered trademark of CytoCell Ltd.

This product contains technology licensed from Life Technologies Corporation that is available for human diagnostics or life science research use only.

### CytoCell Ltd.

Oxford Gene Technology,

418 Cambridge Science Park,

Milton Road,

Cambridge, CB4 0PZ, UK

T: +44(0)1223 294048

F: +44(0)1223 294986

E: probes@cytozell.com

W: www.ogt.com

