



A Sysmex Group Company



TCRAD Breakpart Probe


REF: LPH 047-A / LPH 047-A50

Analyte Specific Reagent: Analytical and performance characteristics are not established.

Fluorescence *in situ* hybridization (FISH) is a technique that allows the visualization of DNA sequences upon chromosomes. The technique uses DNA probes that hybridize to entire chromosomes or single unique sequences, and serves as a powerful adjunct to classic cytogenetics. Recent developments have meant that this valuable technique can now be applied as an essential tool in prenatal, hematological and pathological chromosomal analysis. Target DNA, after fixation and denaturation, is available for annealing to a similarly denatured, fluorescently labelled DNA probe, which has a complementary sequence. Following hybridization, unbound and non-specifically bound DNA probe is removed and the DNA is counterstained for visualization. Fluorescence microscopy then allows the visualization of the hybridized probe on the target material.

Probe Specification

TCRAD, 14q11.2, Red
TCRAD, 14q11.2, Green

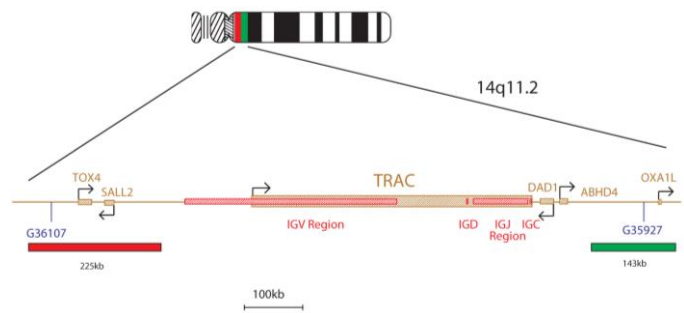

The TCRAD product consists of a 225kb probe, labeled in red, located at the centromeric end of the Variable Region of the immunoglobulin gene cluster, including the G36107 marker and a green probe covering a 143kb region located telomeric to the TRAC gene, including the entire OXA1L gene and the G35927 marker.

Materials Provided

Probe: 100µl per vial or 500µl per vial
Probe concentration: Amount of red probe 1.05-1.77ng/µl
Amount of green probe 6.17-9.23ng/µl

The probe is provided in hybridization solution (Formamide; Dextran Sulphate; SSC) and is ready to use.

Warnings and Precautions

1. For professional use only.
2. Wear gloves when handling DNA probes.
3. Probe contains formamide, which is a teratogen; do not breathe fumes or allow skin contact. Wear gloves, a lab coat, and handle in a fume hood. Upon disposal, flush with a large volume of water.
4. Dispose of all hazardous materials according to your institution's guidelines for hazardous waste disposal.
5. Operators must be capable of visually distinguishing between red, blue and green.

Storage and Handling

The kit should be stored between -25°C to -15°C in a freezer until the expiry date indicated on the kit label. Store the probe vial in the dark. Ensure that exposure of the probe to laboratory lights is limited at all times.

Known Cross-Reactivity

No known cross-reactivity.

Additional Information

For additional product information please contact the CytoCell Technical Support Department.

T: +44 (0)1223 294048
E: techsupport@cytoCELL.com
W: www.ogt.com

Labelling according to Regulation (EC) No. 1272/2008 [CLP]

Hazard pictograms (CLP):



GHS07

GHS08

Signal word (CLP): Danger

Hazardous ingredients: Formamide <100%

Hazard statements (CLP):

H315 – Causes skin irritation
H319 – Causes serious eye irritation
H360 – May damage fertility or the unborn child

Precautionary statements (CLP):

P202 – Do not handle until all safety precautions have been read and understood
P280 – Wear eye protection, protective clothing, protective gloves
P302+P352 – IF ON SKIN: Wash with plenty of soap and water
P305+P351+P338 – IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing
P308+P313 – IF exposed or concerned: Get medical advice/attention
P362+P364 – Take off contaminated clothing and wash it before reuse
P501 – Dispose of contents/container to hazardous or special waste collection point, in accordance with local, regional, national and/or international regulation

Refer to Safety Data Sheet for more information.

Patents and Trademarks

CytoCell is a registered trademark of CytoCell Ltd.
This product contains technology licensed from Life Technologies Corporation that is available for human diagnostics or life science research use only.

CytoCell Ltd.

Oxford Gene Technology,
418 Cambridge Science Park,
Milton Road,
Cambridge, CB4 0PZ, UK
T: +44(0)1223 294048
F: +44(0)1223 294986
E: probes@cytoCELL.com
W: www.ogt.com

