



A Sysmex Group Company



### Instructions For Use

REF: LPH 044-S / LPH 044

## IGH/MAFB Translocation, Dual Fusion Probe



PROFESSIONAL USE ONLY



www.cytoCELL.com

Further information and other languages available at [www.ogt.com](http://www.ogt.com)

### Limitations

This device is designed to detect rearrangements with breakpoints in the region bounded by the red and green clones in this probe set, which includes the *IGH* and *MAFB* genes. Breakpoints outside of this region, or variant rearrangements wholly contained within this region, may not be detected with this product.

The test is not intended for: use as a stand-alone diagnostic, prenatal testing, population-based screening, near-patient testing or self-testing. This product is intended for laboratory professional use only; all results should be interpreted by suitably-qualified staff, taking into account other relevant test results.

This product has not been validated for use on sample types or disease types other than those specified in the intended use.

Reporting and interpretation of FISH results should be consistent with professional standards of practice and should take into consideration other clinical and diagnostic information. This kit is intended as an adjunct to other diagnostic laboratory tests and therapeutic action should not be initiated on the basis of the FISH result alone.

Failure to adhere to the protocol may affect the performance and lead to false positive/negative results.

This kit has not been validated for purposes outside of the intended use stated.

### Intended Use

The CytoCell IGH/MAFB Translocation, Dual Fusion Probe is a qualitative, non-automated, fluorescence *in situ* hybridisation (FISH) test used to detect chromosomal rearrangements between the 14q32.3 region on chromosome 14 and the 20q12 region on chromosome 20 in Carnoy's solution (3:1 methanol/acetic acid) fixed haematologically-derived cell suspensions from patients with confirmed or suspected haematopoietic neoplasias of lymphoid origin [e.g. multiple myeloma (MM)].

### Indications

This product is designed as an adjunct to other clinical and histopathological tests in recognised diagnostic and clinical care pathways, where knowledge of *IGH::MAFB* translocation status would be important for clinical management.

### Principles of the test

Fluorescence *in situ* hybridisation (FISH) is a technique that allows DNA sequences to be detected on metaphase chromosomes or in interphase nuclei from fixed cytogenetic samples. The technique uses DNA probes that hybridise to entire chromosomes or single unique sequences, and serves as a powerful adjunct to G-banded cytogenetic analysis. This technique can now be applied as an essential investigative tool within prenatal, haematological and solid tumour chromosomal analysis. Target DNA, after fixation and denaturation, is available for annealing to a similarly denatured, fluorescently labelled DNA probe, which has a complementary sequence. Following hybridisation, unbound and non-specifically bound DNA probe is removed and the DNA is counterstained for visualisation. Fluorescence microscopy then allows the visualisation of the hybridised probe on the target material.

### Probe Information

The *MAFB* (*MAF bZIP transcription factor B*) gene is located at 20q12 and *IGH* (*immunoglobulin heavy locus*) at 14q32.3

Approximately 50-60% of multiple myeloma (MM) cases are associated with translocations involving *IGH* and one of several partners including *CCND1*, *NSD2* (*MMSET*) and *FGFR3*, *CCND3*, *MAF* or *MAFB*<sup>1</sup>.

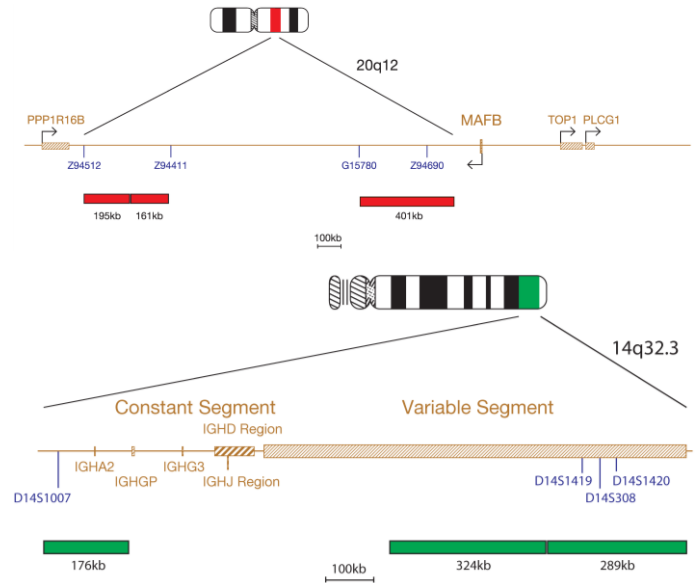
The t(14;20)(q32.3;q12) translocation is a recurrent translocation seen in around 2% of MMs<sup>2,3</sup>.

The reciprocal rearrangement brings a truncated form of the *IGH*  $\mu$ -enhancer ( $E_{\mu}$ , located between the joining (J) segments and the constant region of the *IGH* gene) in close contact with the *MAFB* gene<sup>4</sup>. The resultant fusion and the up-regulated transcription product has been shown to cause dysregulation of cyclin D2<sup>1</sup>.

The prognostic outcome of t(14;20)(q32.3;q12) is assumed to be the same as the t(14;16)(q32.3;q23)<sup>3</sup>.

### Probe Specification

MAFB, 20q12, Red  
IGH, 14q32.3, Green



The IGH/MAFB product consists of probes, labelled in green, covering the Constant and Variable segments of the *IGH* gene, and *MAFB* probes (195kb, 161kb and 401kb), labelled in red. The *MAFB* probes are located on either side of the breakpoint region (between *MAFB* and *PPP1R16B*).

### Materials Provided

**Probe:** 50 $\mu$ l per vial (5 tests) or 100 $\mu$ l per vial (10 tests)

The probes are provided premixed in hybridisation solution (formamide; dextran sulphate; saline-sodium citrate (SSC)) and are ready to use.

**Counterstain:** 150 $\mu$ l per vial (15 tests)

The counterstain is DAPI antifade (ES: 0.125 $\mu$ g/ml DAPI (4,6-diamidino-2-phenylindole)).

### Warnings and Precautions

1. For *in vitro* diagnostic use. For professional use only.
2. Wear gloves when handling DNA probes and DAPI counterstain.
3. Probe mixtures contain formamide, which is a teratogen; do not breathe fumes or allow skin contact. Wear gloves, a lab coat, and handle in a fume hood. Upon disposal, flush with a large volume of water.
4. DAPI is a potential carcinogen. Handle with care; wear gloves and a lab coat. Upon disposal, flush with a large volume of water.
5. Dispose of all hazardous materials according to your institution's guidelines for hazardous waste disposal.
6. Operators must be capable of distinguishing the colours red, blue and green.
7. Failure to adhere to the outlined protocol and reagents may affect the performance and lead to false positive/negative results.
8. The probe should not be diluted or mixed with other probes.
9. Failure to use 10 $\mu$ l of probe during the pre-denaturation stage of the protocol may affect the performance and lead to false positive/negative results.

### Storage and Handling

The kit should be stored between -25 $^{\circ}$ C to -15 $^{\circ}$ C in a freezer until the expiry date indicated on the kit label. The probe and counterstain vials must be stored in the dark.



The probe remains stable throughout the freeze-thaw cycles experienced during normal use (where one cycle constitutes the probe's removal from and replacement into the freezer) and is photostable for up to 48 hours after being exposed to continuous lighting conditions. All efforts must be made to limit exposure to light and temperature changes.

## Equipment and Materials Necessary but not Supplied

Calibrated equipment must be used:

1. Hotplate (with a solid plate and accurate temperature control up to 80°C)
2. Calibrated variable volume micropipettes and tips range 1µl - 200µl
3. Water bath with accurate temperature control at 37°C and 72°C
4. Microcentrifuge tubes (0.5ml)
5. Fluorescence microscope (Please see Fluorescence Microscope Recommendation section)
6. Phase contrast microscope
7. Clean plastic, ceramic or heat-resistant glass Coplin jars
8. Forceps
9. Calibrated pH meter (or pH indicator strips capable of measuring pH 6.5 – 8.0)
10. Humidified container
11. Fluorescence grade microscope lens immersion oil
12. Bench top centrifuge
13. Microscope slides
14. 24x24mm coverslips
15. Timer
16. 37°C incubator
17. Rubber solution glue
18. Vortex mixer
19. Graduated cylinders
20. Magnetic stirrer
21. Calibrated thermometer

## Optional Equipment not Supplied

1. Cytogenetic drying chamber

## Reagents Needed but not Supplied

1. 20x saline-sodium citrate (SSC) Solution
2. 100% Ethanol
3. Tween-20
4. 1M Sodium hydroxide (NaOH)
5. 1M Hydrochloric acid (HCl)
6. Purified water

## Fluorescence Microscope Recommendation

Use a 100-watt mercury lamp or equivalent and oil immersion plan apochromat objectives 60/63x or 100x for optimal visualisation. The fluorophores used in this probe set will excite and emit at the following wavelengths:

| Fluorophore | Excitation <sub>max</sub> [nm] | Emission <sub>max</sub> [nm] |
|-------------|--------------------------------|------------------------------|
| Green       | 495                            | 521                          |
| Red         | 596                            | 615                          |

Ensure appropriate excitation and emission filters that cover the wavelengths listed above are fitted to the microscope. Use a triple bandpass DAPI/green spectrum/red spectrum filter or a dual bandpass green spectrum/red spectrum filter for optimal simultaneous visualisation of the green and red fluorophores.

Check the fluorescence microscope before use to ensure it is operating correctly. Use immersion oil that is suitable for fluorescence microscopy and formulated for low auto fluorescence. Avoid mixing DAPI antifade with microscope immersion oil as this will obscure signals. Follow manufacturers' recommendations in regards to the life of the lamp and the age of the filters.

## Sample Preparation

The kit is designed for use on haematologically-derived cell suspensions fixed in Carnoy's solution (3:1 methanol/acetic acid) fixative, that are prepared according to the laboratory or institution guidelines. Prepare air dried samples on microscope slides according to standard cytogenetic procedures. The AGT *Cytogenetics Laboratory Manual* contains recommendations for specimen collection, culturing, harvesting and for slide making<sup>5</sup>.

## Solution Preparation

### Ethanol Solutions

Dilute 100% ethanol with purified water using the following ratios and mix thoroughly:

- 70% Ethanol - 7 parts 100% ethanol to 3 parts purified water
- 85% Ethanol - 8.5 parts 100% ethanol to 1.5 parts purified water

Store the solutions for up to 6 months at room temperature in an airtight container.

### 2xSSC Solution

Dilute 1 part 20xSSC Solution with 9 parts purified water and mix thoroughly. Check pH and adjust to pH 7.0 using NaOH or HCl as required. Store the solution for up to 4 weeks at room temperature in an airtight container.

### 0.4xSSC Solution

Dilute 1 part 20xSSC Solution with 49 parts purified water and mix thoroughly. Check pH and adjust to pH 7.0 using NaOH or HCl as required. Store the solution for up to 4 weeks at room temperature in an airtight container.

### 2xSSC, 0.05% Tween-20 Solution

Dilute 1 part 20xSSC Solution with 9 parts purified water. Add 5µl of Tween-20 per 10ml and mix thoroughly. Check pH and adjust to pH 7.0 using NaOH or HCl as required. Store the solution for up to 4 weeks at room temperature in an airtight container.

## FISH Protocol

(Note: Ensure that exposure of the probe and counterstain to laboratory lights is limited at all times).

### Slide Preparation

1. Spot the cell sample onto a glass microscope slide. Allow to dry. (**Optional, if using a cytogenetic drying chamber**: slides should be spotted using a cytogenetic drying chamber. The chamber should be operated at approximately 25°C and 50% humidity for optimal cell sample spotting. If a cytogenetic drying chamber is not available, use a fume hood as an alternative).
2. Immerse the slide in 2xSSC for 2 minutes at room temperature (RT) without agitation.
3. Dehydrate in an ethanol series (70%, 85% and 100%), each for 2 minutes at RT.
4. Allow to dry.

### Pre-Denaturation

5. Remove the probe from the freezer and allow it to warm to RT. Briefly centrifuge tubes before use.
6. Ensure that the probe solution is uniformly mixed with a pipette.
7. Remove 10µl of probe per test, and transfer it to a microcentrifuge tube. Quickly return the remaining probe to the freezer.
8. Place the probe and the sample slide to prewarm on a 37°C (+/- 1°C) hotplate for 5 minutes.
9. Spot 10µl of probe mixture onto the cell sample and carefully apply a coverslip. Seal with rubber solution glue and allow the glue to dry completely.

### Denaturation

10. Denature the sample and probe simultaneously by heating the slide on a hotplate at 75°C (+/- 1°C) for 2 minutes.

### Hybridisation

11. Place the slide in a humid, lightproof container at 37°C (+/- 1°C) overnight.

### Post-Hybridisation Washes

12. Remove the DAPI from the freezer and allow it to warm to RT.
13. Remove the coverslip and all traces of glue carefully.
14. Immerse the slide in 0.4xSSC (pH 7.0) at 72°C (+/- 1°C) for 2 minutes without agitation.
15. Drain the slide and immerse it in 2xSSC, 0.05% Tween-20 at RT (pH 7.0) for 30 seconds without agitation.
16. Drain the slide and apply 10µl of DAPI antifade onto each sample.
17. Cover with a coverslip, remove any bubbles and allow the colour to develop in the dark for 10 minutes.
18. View with a fluorescence microscope (see **Fluorescence Microscope Recommendation**).

### Stability of Finished Slides

Finished slides remain analysable for up to 1 month if stored in the dark at/or below RT.

### Procedural Recommendations

1. Baking or ageing of slides may reduce signal fluorescence
2. Hybridisation conditions may be adversely affected by the use of reagents other than those provided or recommended by Cytocell Ltd
3. Use a calibrated thermometer for measuring temperatures of solutions, waterbaths and incubators as these temperatures are critical for optimum product performance.
4. The wash concentrations, pH and temperatures are important as low stringency can result in non-specific binding of the probe and too high stringency can result in a lack of signal
5. Incomplete denaturation can result in lack of signal and over denaturation can also result in non-specific binding
6. Over hybridisation can result in additional or unexpected signals
7. Users should optimise the protocol for their own samples prior to using the test for diagnostic purposes
8. Suboptimal conditions may result in non-specific binding that may be misinterpreted as a probe signal

### Interpretation of Results

#### Assessing Slide Quality

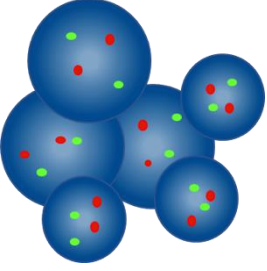
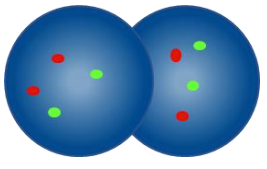
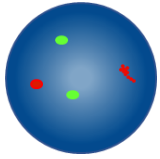
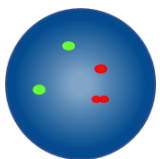
The slide should not be analysed if:

- Signals are too weak to analyse in single filters - in order to proceed with analysis, signals should appear bright, distinct and easily evaluable
- There are high numbers of clumped/overlapping cells obstructing the analysis
- >50% of the cells are not hybridised
- There is excess of fluorescent particles between cells and/or a fluorescent haze that interferes with the signals - in optimal slides the background should appear dark or black and clean
- Cell nucleus borders cannot be distinguished and are not intact

### Analysis Guidelines


- Two analysts should analyse and interpret each sample. Any discrepancies should be resolved by assessment by a third analyst
- Each analyst should be suitably qualified according to recognised national standards
- Each analyst should score independently 100 nuclei for each sample. The first analyst should start the analysis from the left side of the slide and the second analyst from the right one
- Each analyst should document their results in separate sheets

- Analyse only intact nuclei, not overlapped or crowded nuclei or nuclei covered by cytoplasmic debris or high degree of autofluorescence
- Avoid areas where there is excess of cytoplasmic debris or non-specific hybridisation
- Signal intensity may vary, even with a single nucleus. In such cases, use single filters and/or adjust the focal plane
- In suboptimal conditions signals may appear diffuse. If two signals of the same colour touch each other, or the distance between them is no greater than two signal widths, or when there is a faint strand connecting the two signals, count as one signal
- If in doubt about whether a cell is analysable or not, then do not analyse it

| Analysis Guidelines   |   |
|---|---|
|    | Do not count – nuclei are too close together to determine boundaries                                      |
|    | Do not count overlapping nuclei – all areas of both nuclei are not visible                                |
|   | Count as two red signals and two green signals – one of the two red signals is diffuse                    |
|  | Count as two red signals and two green signals – the gap in one red signal is less than two signal widths |


#### Expected Results
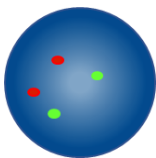
##### Expected normal signal pattern



In a normal cell, two red and two green signals (2R, 2G) are expected.

##### Expected abnormal signal pattern

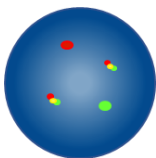

In a cell with a t(14;20)(q32.3;q12) translocation, one red, one green and two red/green fusion signals (1R, 1G, 2F) are expected.

Other signal patterns are possible in aneuploid/unbalanced specimens.

#### Known Cross-Reactivity

The green IGH probe may show cross-hybridisation to 15q11.2 and 16p11.2.

#### Adverse Event Reporting

If you believe this device has malfunctioned or suffered a deterioration in its performance characteristics which may have contributed to an adverse event (e.g. delayed or misdiagnosis, delayed or inappropriate treatment), this must be reported immediately to the manufacturer (**email:** [vigilance@ogt.com](mailto:vigilance@ogt.com)).

If applicable, the event should also be reported to your national competent authority. A list of vigilance contact points can be found at: <http://ec.europa.eu/growth/sectors/medical-devices/contacts/>.

#### Specific Performance Characteristics

Laboratories must verify cut-off values using their own data<sup>6,7</sup>.

#### Analytical Specificity

Analytical specificity is the percentage of signals that hybridise to the correct locus and no other location. The analytical specificity was established by analysing a total of 200 target loci. The analytical specificity was calculated as the number of FISH signals that hybridised to the correct locus divided by the total number of FISH signals hybridised.

Table 1. Analytical Specificity for the IGH/MAFB Translocation, Dual Fusion Probe

| Probe     | Target Locus | No. of Signals Hybridised to the Correct Locus | Total No. of Signals Hybridised | Specificity (%) |
|-----------|--------------|--|---------------------------------|-----------------|
| Red MAFB  | 20q12        | 200  | 200                             | 100             |
| Green IGH | 14q32.33     | 200  | 200                             | 100             |

#### Analytical Sensitivity

Analytical sensitivity is the percentage of scoreable interphase cells with the expected normal signal pattern. The analytical sensitivity was established by analysing interphase cells across different normal samples. The sensitivity was calculated as the percentage of scoreable cells with the expected signal pattern (with a 95% confidence interval).

Table 2. Analytical Sensitivity for the IGH/MAFB Translocation, Dual Fusion Probe

| No. of Cells with Expected Signal Patterns | No. of Cells with Scoreable Signals | Sensitivity (%) | 95% Confidence Interval |
|--|-------------------------------------|-----------------|-------------------------|
| 467  | 500                                 | 93.4            | 2.1                     |

#### Precision and Reproducibility

Precision is a measure of the natural variation of a test when repeated several times under the same conditions. This was assessed by analysing repeats of the same lot number of probe tested on the same sample, in the same conditions on the same day.

Reproducibility is a measure of the variability of a test and has been established in terms of sample-to-sample, day-to-day and batch-to-batch variability. Day-to-day reproducibility was assessed by analysing the same samples on three different days. Batch-to-batch reproducibility was assessed by analysing the same samples using three different lot numbers of probe on one day. Sample-to-sample reproducibility was assessed by analysing three replicates of a sample on one day. For each sample, signal patterns of 100 interphase cells were recorded and the percentage of cells with the expected signal pattern was calculated.

The reproducibility and precision were calculated as the Standard Deviation (STDEV) between replicates for each variable and overall mean STDEV.

Table 3. Reproducibility and Precision for the IGH/MAFB Translocation, Dual Fusion Probe

| Variable          | Standard Deviation (STDEV) |
|-------------------|----------------------------|
| Precision         | 0.00                       |
| Sample-to-sample  | 0.00                       |
| Day-to-day        | 0.00                       |
| Batch-to-batch    | 0.00                       |
| Overall deviation | 0.00                       |

#### Additional Information

For additional product information please contact the CytoCell Technical Support Department.

T: +44 (0)1223 294048


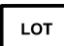



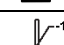



E: [techsupport@cytoCELL.com](mailto:techsupport@cytoCELL.com)

W: [www.ogt.com](http://www.ogt.com)

#### References

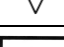- Fonseca *et al.*, Cancer Research 2004;64:1546-1558
- Fonseca *et al.*, Leukemia 2009;23(12):2210-2221
- Sawyer, Cancer Genetics 2011;204(1):3-12
- Boersma-Vreugdenhil *et al.*, Br J Haematol 2004;126:355-63
- Arsham, MS., Barch, M.J. and Lawce H.J. (eds.) (2017) *The AGT Cytogenetics Laboratory Manual*. New Jersey: John Wiley & Sons Inc.
- Mascarello JT, Hirsch B, Kearney HM, et al. Section E9 of the American College of Medical Genetics technical standards and guidelines: fluorescence in situ hybridization. Genet Med. 2011;13(7):667-675.
- Wiktor AE, Dyke DLV, Stupca PJ, Ketterling RP, Thorland EC, Shearer BM, Fink SR, Stockero KJ, Majorowicz JR, Dewald GW. *Preclinical validation of fluorescence in situ hybridization assays for clinical practice*. Genetics in Medicine. 2006;8(1):16–23.

## Guide to Symbols

|   |   |
|---|---|
| REF   | en: Catalogue number                          |
|  | en: <i>In vitro</i> diagnostic medical device |
|  | en: Batch code                                |
|  | en: Consult instructions for use              |
|  | en: Manufacturer                              |
|  | en: Use-by date                               |
|  | en: Temperature limit                         |
|  | en: Keep away from sunlight                   |
|  | en: Contains sufficient for <n> tests         |
|  | en: Contents                                  |

## Patents and Trademarks

CytoCell is a registered trademark of CytoCell Ltd.



### Cytocell Ltd.

Oxford Gene Technology,  
418 Cambridge Science Park,  
Milton Road,  
Cambridge, CB4 0PZ, UK  
T: +44(0)1223 294048  
F: +44(0)1223 294986  
E: probes@cytoCell.com  
W: www.ogt.com